



Customer							
Material	Placebo Beads 0.8mm dia.		Est. Bulk Density	3.0 gm/cc EST.		Date	4/26/02
Desired Sample size	Various capsules sizes #0 to #4 1/4 to 3/4 fill						
Desired accuracy	NA		Pipette Size	.250 dia.		rh	50%

Micrometer setting	Sample weight	Control unit settings	
	mg	Vacuum (1in hg)	Air (psi)
1.50	0.0195	5	2-Jan
#4 capsule	0.0201		
	0.0261		
	0.0232		
	0.0185		
	0.0202		
	0.0244		
	0.0231		
	0.0225		
	0.0280		
	0.0223		
		0.0225 Av	
	0.0240 Mean		
	0.0029 Std Dev		

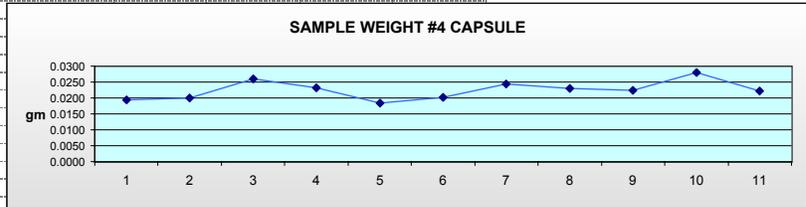

4.00	0.2611			
#2 CAPSULE	0.2624			
	0.2567			
	0.2579			
	0.2622			
	0.2624			
	0.2602			
	0.2607			
	0.2610			
	0.2614			
		0.2606 Av		
		0.2606 Mean		
	0.0019 Std Dev			

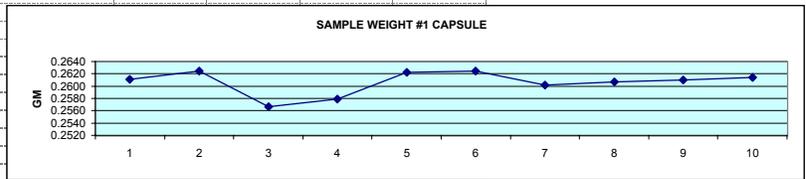

9.0	0.5595			
#0 CAPSULE	0.5587			
	0.5560			
	0.5569			
	0.5602			
	0.5612			
	0.5591			
	0.5635			
	0.5633			
	0.5629			
		0.5601 Av		
		0.5601 Mean		
	0.0026 Std Dev			

